

Invented for life



**BOSCH**









# Technical Datasheet

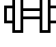





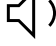

√ivalytic from Bosch –  
The universal diagnostic platform.  
**Universal, intuitive, fast.**




√ivalytic

# Technical Datasheet

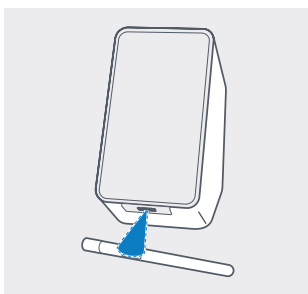
Description	Value
 Display	7 inch 16:10, 1024 x 600 pixel touchscreen
 Operating air pressure range	850-1,100 hPa, corresponds to pressure range 0-1, 400 m above sea level
 Operating temperature	15-30 °C
 Storage temperature	-20-60 °C
 Data transfer	Ethernet 10/100 MB, WLAN 2.4 GHz (802.11b/g/n); internal: Bluetooth v4.1, 2.4 GHz (low energy), USB 2.0
 Electromagnetic compatibility	IEC/EN 61326-2-6 RED 2014/53/EC FCC47 CFR 15
 Dimensions	Length 400 mm Width 204 mm Height 388 mm
 Distance/Space to the wall	~20 cm

Description	Value
 Weight	15 kg
 Storage humidity	20-95 % (not condensing)
 Operating humidity	30-80 % (not condensing)
 Electrical data	100-240 V~, 50/60 Hz, 160 VA
 Instrumental safety	IEC/EN 61010-1 IEC/EN 61010-2-010 IEC/EN 61010-2-101 Regulation 2017/746
 Memory capacity	16 GB
 Mean loudness	≤ 55 dB(A) in operating mode. Short term loudness can exceed mean loudness
 Socket	Use multiple sockets for EU countries and UK


## Vivalytic Workflow: Good results can be this easy.




1 Scan sample code.

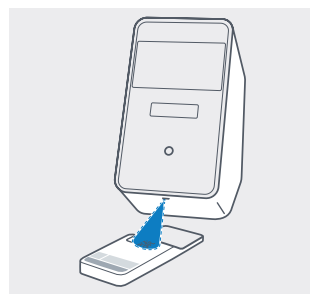

2 Scan cartridge code.

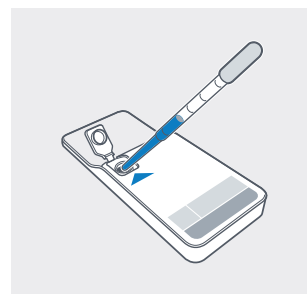

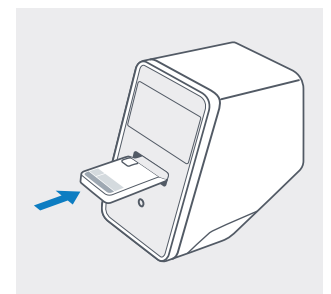3 Insert the sample and close cartridge lid.



4 Insert cartridge. The test starts automatically.

Subject to technical modifications. Printed in Germany. Any colour deviations are due to restrictions in the printing process. Printed on 100% chlorine-free bleached cellulose, produced in an environmentally friendly process. Not all products are available in all regions. Ask your local sales representative for availability in specific markets. For in vitro diagnostic use only.

**Bosch Healthcare Solutions GmbH**  
Stuttgarter Str. 130, 71332 Waiblingen, Germany